

Welcome to Foot Pedal Prediction!

Press any button to continue.

Please take note of the control
around your right foot.

Hit the up button to continue.

Please take note of the control
around your right foot.

Hit the right button to continue.


Please take note of the
control around your left foot.

Hit the left button to continue.

Great job! We are an application
designed to help Brad
communicate more seamlessly
with those around him.

Hit the any button to continue.

When you navigate to the main application page, you will see something like this:

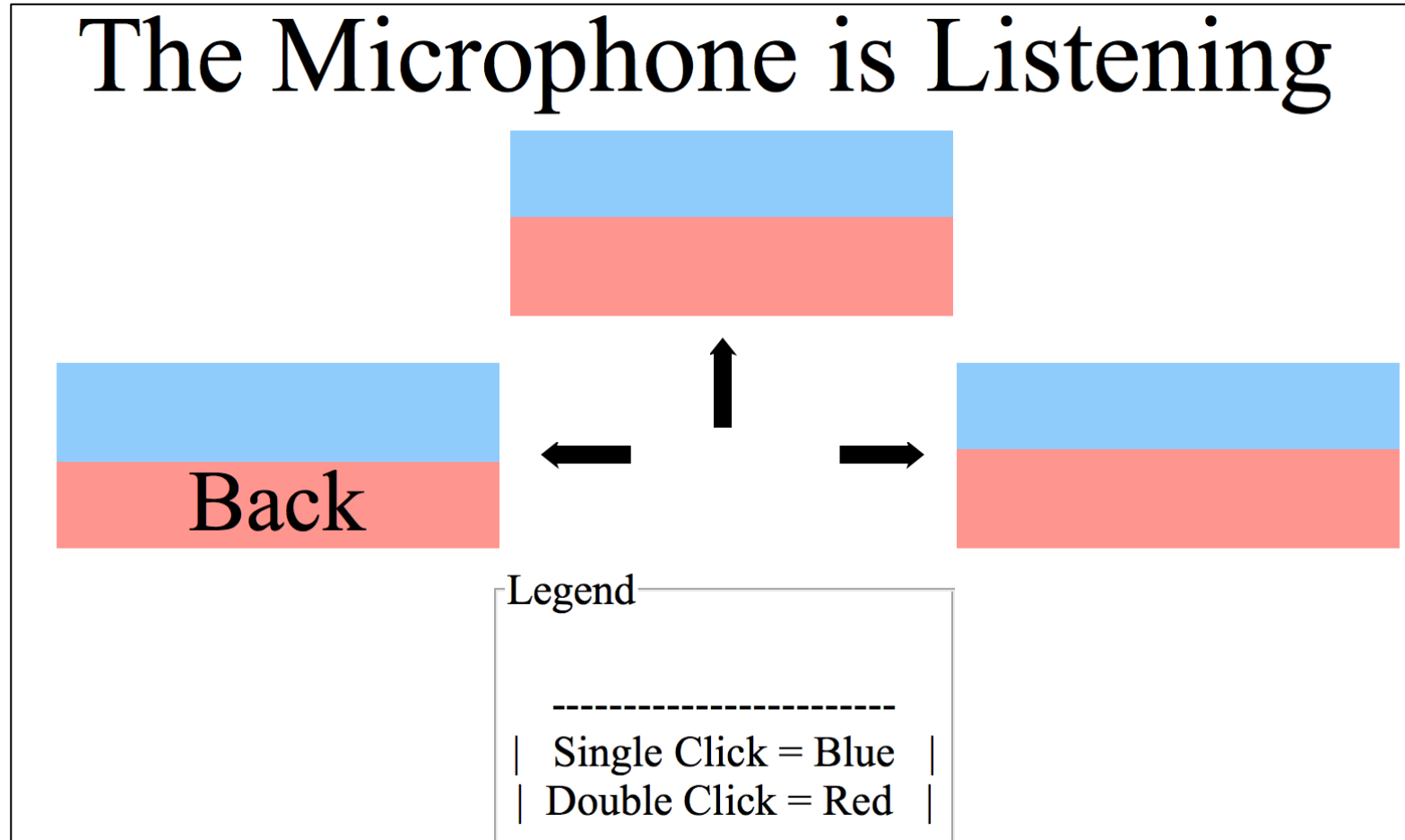

Hit any button to continue.

As you speak, the colored squares will fill with words we predict you may say next. You can select the **blue words** with a single click, and the **red words** with a double click. Let's practice a double click now.

Double-click any button to continue.

The **red** button labelled “Back” will always be positioned in its place. A double-click on the left button from the main application will return you to the main menu. Let’s practice now.

Double-click the left button to continue.

When you navigate to the settings page,
you will see something like this:

Settings

Press the up key to rotate.
Left to lessen variable.
Right to increase variable value.

Number of words for the app to display is: 3

Number of seconds for the microphone to sleep after listening: 3

Number of seconds to consider the double-click interval: 1

Exit

Hit any button to continue.

When on the settings page, use the up bottom to rotate through the selected setting, and the left and right arrow keys to adjust them according to your preference (pressing either with “Exit” highlighted will return you to the main menu).

Click any button to continue.

Thank you and we hope you
enjoy using our product!

Press any button to return to the main menu.